



# Molift QuickRaiser 205

A sit to stand hoist with outstanding  
lifting capacity

# Molift QuickRaiser 205

A sit to stand hoist with outstanding lifting capacity

Molift QuickRaiser 205 is a state of the art sit-to-stand hoist, with an electrically adjustable leg base and an outstanding lifting capacity. The standing hoist is comfortable, stable and safe with exceptional manoeuvrability.

## Stable and compact

QuickRaiser 205 is compact, stable and exceptionally easy to manoeuvre. The unique wheel position, angled push bar and the extraordinary weight distribution of the base provides a small turning radius and easy manoeuvrability even in very narrow spaces. With the small castor option (50 mm), it is only 7,5 cm (3") high from the floor, and fits under most beds and chairs while the electrically adjustable legs enables accessibility almost everywhere.

## Natural pattern of movement

The Active Lifting Arm (four point spreader bar) and the RgoSling Active provides extra support when hoisting from sit to stand. The inclined column encourages a natural pattern of movement and fulfills the hoisting needs of users with some weight bearing capacity. The Active Lifting Arm combined with the sling provides a comfortable, stable and safe hoisting solution.



Molift QuickRaiser 205

Description	Castor	Item no.	Width of wheelbase	
Molift QuickRaiser 205 Active	50/75mm	29000	External	62 cm
Molift QuickRaiser 205 Standup	50/75mm	29001	Internal	20 cm
Molift QuickRaiser 205 Active	75/100mm	29100		
Molift QuickRaiser 205 Standup	75/100mm	29101		
Weight			Castor	Height of wheelbase
Maximum user weight (SWL)		205 kg	Ø 75/50mm	76 mm
Total weight		44 kg	Ø 100/75mm	101 mm
Chassis		17 kg		
Lifting arm		4,5 kg		
Leg support		3,5 kg		
Column (heaviest part)		19 kg		
Battery SLA 24V		2,9 Ah	Castor	Ground clearance
Charge time		6 hours	Ø 75/50mm	27 mm
Number of lifts (75 kg) with fully charged battery		>100	Ø 100/75mm	52 mm
Speed			Castor	Height of footplate
Hoisting speed 75 kg		3,5 cm/s	Ø 75/50mm	65 mm
Hoisting speed 205 kg		3,0 cm/s	Ø 100/75mm	90 mm
Lifting Range		60 cm	Turning radius	115 cm
			Cleaning instructions	
			Clean the product with a solvent-free cleaning agent (pH level 5–9), or with a 70% disinfectant solution. Can be decontaminated at max 85°C / 185°F	



## Features



### Ergonomic design

The slanting column lifts the user up and forwards to imitate a natural movement pattern.



### Automatic Warning

Automatic warning for (annual) periodic inspection and service.



### Comfort, safety and stability

The Active Lifting Arm (four point spreader bar) and the RgoSling Active provides extra support when hoisting from sit to stand.




### Easy to manoeuver

Easy to operate in moving and handling situations; toilet, sit-to-stand and stand-up training.



*The new operation panel and hand control provides intuitive and safe operation. Furthermore it includes*

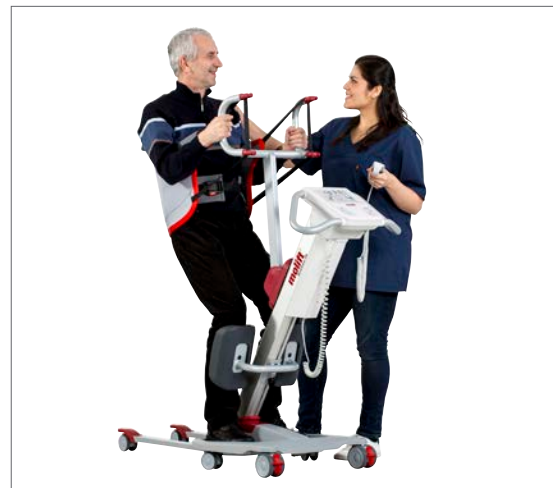

*Active lifting Arm and RgoSling Active.*

## Accessories



### StandUp lifting arm

Ideal to use with RgoSling Standup slings.

StandUp Lifting Arm 2910061



### Inlays for footplate

Inlays with slipresistant material to increase friction.

Set of inlays for footplate 2920151



*"Everyone should be able to live a free and independent life  
and pursue their dreams regardless of any physical circumstances."*

*For the latest news and continuously  
updated product information – please  
visit **www.etac.com***

Etac Supply Gjøvik  
Etac AS  
Hadelandsveien 2  
N-2816 Gjøvik  
molift@etac.com

---

Etac AB  
Kista Science Tower  
SE-164 51 Kista, Sweden  
Tel +46 371 – 58 73 30  
Fax +46 371 – 58 73 90  
info@etac.se www.etac.com

**etac**<sup>®</sup>  
Creating Possibilities