

Etac Swift Mobile Tilt

With a unique tilt in space function.

Etac Swift Mobile Tilt is available in two models. Both have manually operated tilt functions but with different weight capacity. They all have the same characteristic design and functions as Etac Swift Mobile, page 26. As an option Swift Mobile Tilt can be equipped with a power unit for the tilt function.

Tilting operated from the side

Etac Swift Mobile Tilt has the unique function which enables tilt manoeuvring from the side. This gives the best conditions for an ergonomically correct work position for the carer. It also enables eye contact and communication with the care recipient during the tilting process.

Assisted restoring to an upright position

The gas pistons work to facilitate bringing the seat into an upright position. When tilting, the user's body weight will act as counterweight and make the movement smooth. This reduces the risk of carers suffering work related injuries. The seat height can be set in various fixed positions and it's done without using tools.



A tilted position increases stability and facilitates in care situations.



Power tilt function for easy and smooth manoeuvre

Etac Swift Mobile Tilt can be equipped with a power unit to control the tilt function which is maneuvered with a practical hand control. This means that the user is also able to control the tilt function. The Power unit occupy very little space, which means that the shower stool can still be positioned over the toilet. The chair is approved to be used in the shower and a separate battery charger is included.

Features



Unique tilt bar

The bar for the tilt mechanism is easy to reach from either side.



Head Support

The adjustability ensures a comfortable and restful support.





Close contact during the tilting process

The unique tilt function enables an ergonomically correct working position for the carer who is able to maintain eye contact with the care recipient.



Power tilt unit

Enables a better working environment. Smooth manoeuvre with a practical hand control.

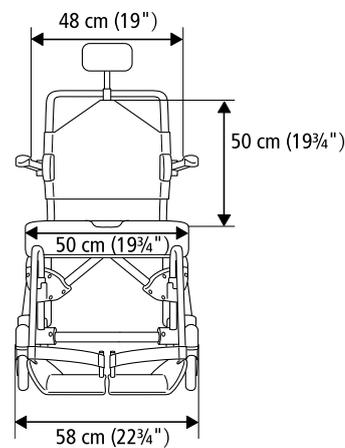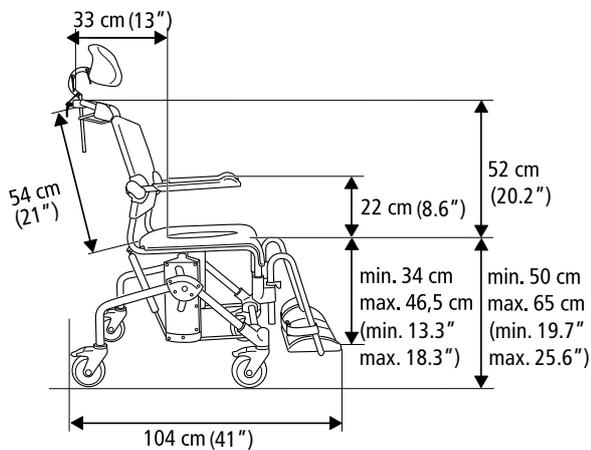

Accessories

See our wide range of accessories
"Accessories for all Swift Mobile models"
on page 34

Description Etac Swift Mobile Tilt

Swift Mobile Tilt	80209410
Swift Mobile Tilt with bucket holder	80209411
Swift Mobile Tilt with pan holder	80209412
Swift Mobile Tilt with power unit	80209490
Swift Mobile Tilt with power unit incl. bucket holder	80209491
Swift Mobile Tilt with power unit incl. pan holder	80209492

Standard: head support, back support, arm supports, seat, foot supports, heels traps, and 4 lockable castors.

Dimensions

Adjustable tilt in space: 0° to -35° can be set to +5° to -30°

Material

Gas pistons: lift capacity 130N
Other material: see "Etac Swift Mobile shower commode"
on page 26

Weight

19.7 kg (43.4 lbs)
With power unit: 23.6 kg (52 lbs)



Clean the product with a solvent-free cleaning agent (pH level 6-9), or with a 70% disinfectant solution. Can be decontaminated at max. 85°C / 185°F if so required.

Accessories for all Swift Mobile models



Head Support

80209447

Complete with attachment. Adjustable in height, depth and sideways. Standard on tilt models.

Material: polyurethane, stainless steel and zinc. Cleaning max. 85° C (185°F) or with alcohol wipe.



Comfort Cover

80209455

A soft back support with integrated arm support pads.

Material: polyethylene.

Cleaning max. 60° C (140 °F) or with alcohol wipe.



Arm supports widening kit

Pair: 80209441

Increases the space between the arm supports by 6 cm (2½"). Standard for Swift Mobile/Tilt 160.



Comfort Kit

80209457

Includes both Comfort Seat and Comfort Cover.

Cleaning max. 60° C (140 °F) or with alcohol wipe.



Hemi Arm Support

Pair: 80209035

Fits both left and right side. Material: Cellular polyethylene, aluminium, stainless steel, polystyrene.



Comfort arm pads – pair

80209508

The pads provide a soft support for the arms.

Material: polyethene

Cleaning max. 60° C (140 °F) or with alcohol wipe.



Arm supports lock

Pair: 80209445

Locks the arm supports in down-folded position. Material: Stainless steel.



Side support

80209509

Adjustable in width.

Material: Polyester fabric, Velcro, polyethylene.

Cleaning at max 60°C (140°F)



Soft Back Pad

80209448

Warm, soft and water repellent.

Material: Polyethylene

Cleaning max. 60° C (140°F) or with alcohol wipe



Soft cushion with standard opening

80209065

Material: Polyether/Sealed coated cover



Seat cover

80209434

When a complete seat is required.

Material: Foamed polyurethane.

max. 85° C (185°F) or with alcohol wipe



Soft cushion with oval opening

80209067

Material: Polyether/Sealed coated cover



Soft Comfort seat

Opening width 15 cm (6"): 80209506

Opening width 18 cm (7"): 80209507

Material: Polyurethane

Cleaning at max 85°C (185°F) or with alcohol wipe.



Chest harness strap

80209442

Adjustable in length. Attach to back support. Always use together with positioning belt.

Material: Polyester



Cushion, Protect

80209245

Air-inflated ROHO cushion with excellent pressure distribution properties. Adjustable air pressure. Recess width: 21 cm (8¼").

Material: Neoprene

Cleaning max: 60° C (140 °F)



Cross bar

Length: 60.5 cm (24"): 80303034

Length: 66.5 cm (26"): 80209452

For safety and additional support.

Material: Aluminium, polypropene and cellular polyethylene. Notice: 66.5 cm is for Swift Mobile 160 and Swift Mobile Tilt 160.



Positioning belt
80209440
Adjustable in length. Attachment with velcro straps.
Material: Polyester



Elevating Leg- and calf support
80209437
Fits both left and right side.
Material: ABS



Splash Guard
80209430
Fits into the recess opening. To be placed when user sits in the chair.
Material: Polyurethane



Foot support for leg- and calf support
80209443
Fits both left and right side. Foot support adjustable height in 6 positions.
Material: Stainless steel, polypropylene.



Attachment for pan/bucket holder
Set: 80209435
To be mounted underneath the seat.
Material: Polyamide.



Leg support
1 piece: 80209446
Leg support with foldable and angle adjustable footplate. No heel strap included.



Pan holder
80209257
Fits both Pan and regular hospital pans with a thin rim. Attachment required.
Material: Polyamide.



Heel straps
Pair: 80209432
Adjustable in length.
Material: Polyester



Pan
80209254



Anti tippers, pair
Excluding attachment for Swift Mobile 24": 80209438
Including attachment (for Swift Mobile): 80209439
(for Swift Mobile Tilt): 80209459
Increases safety if the center of gravity is far back, e.g. leg amputations.
Material: Aluminium



Bed pan
82001001
Polyamide, able to withstand boiling.
Fits the Pan holder.



Tilter
80209449
For easy transfer over thresholds and other obstacles



Bucket holder
80209258
Fits bucket with a thick rim. Attachment required.
Material: Polypropene



Rear wheel kit
80209436
24" rear wheels with solid tires. With built in attachment for anti tip.
Material: Polyamide, stainless steel.
Incl. attachment and brakes. Do not fit Swift Mobile 160 or Swift Mobile Tilt.



Bucket with handle and lid
80302026
Material: Polypropylene



Attachment for low mounting of pan/bucket
80209259
Fitted to the pan/bucket holder.
Material: Polypropene



Power tilt unit
80209499
For Swift Mobile Tilt and Swift Mobile Tilt 160 only. Hand control and charger are included.



Amputee leg support
80209510
Adjustable height, length and angle. Fits both left and right side.