

Proven performance, quality & service. When it counts.

Dear Valued Customer

Our Mission Statement summarizes our corporate philosophy...providing our customers with superior products, education and services which improve quality of life through shape fitting solutions. Our pledge to the customer is to provide innovative products that are durable, lightweight and adjustable, and which ultimately enhance the customer's quality of life.

We are proud of our ROHO products and believe them to be one of the best seating solutions, for every type of person.

We have built our reputation not only on quality products, but also on strong customer service, solid personal relationships, supporting sales and training programs, and an advanced manufacturing process dedicated to continually improving quality for our customer base. We are a company that sees the customer's needs as our first priority.

The ROHO Group greatly appreciates your continued interest and support.
We are committed to providing you with the very best in quality rehabilitation and wound care products.

Sincerely,



TABLE OF CONTENTS

ROHO® DRY FLOATATION® Technology4-6
Seat Cushions
ROHO® QUADTRO SELECT® Cushion 10
ROHO® CONTOUR SELECT® Cushion 11
ROHO® HIGH PROFILE® Cushion12
ROHO® MID PROFILE™ Cushion12
ROHO® LOW PROFILE® Cushion
ROHO® ENHANCER® Cushion13
ROHO Hybrid Elite® Cushion14
ROHO® nexus SPIRIT® Cushion14
ROHO® HARMONY® Cushion15
ROHO® AirLITE® Cushion
ROHO® MOSAIC® Cushion
PostureLITE™ by ROHO® Cushion17
LTV Seat® Cushion
Seat Cushion Sizing Chart18-21
Support Surface Products
ROHO® DRY FLOATATION® Mattress Overlay 26
ROHO® DRY FLOATATION® Bariatric Kit
ROHO® DRY FLOATATION® Isolette
Reusable Mattress Cover27
ROHO® PRODIGY® Mattress Overlay28
ROHO® SOFFLEX® 2 Mattress Overlay System 28
Sierra™ Gel Overlay29
SelectAir® Low Air Loss System30
SelectAir® Max Low Air Loss System31
SelectProtect® Safety Enhancement System 31
Bari Select® Low Air Loss System32
Hybrid <i>Select</i> ® with Heal <i>Float</i> ®
Low Air Loss System33
Fusion 1K [™] Mattress Replacement System 34-35
Fusion 2K [™] Mattress Replacement System 36-37
Fusion XC™ Mattress Replacement System 38-39

Dack Support Systems	
AGILITY™ Minimum Contour Back	42
AGILITY™ Mid Contour Back	43
JetStream Pro® Back Support System	44
Retroback® Back Support System	45
Custom Recliner System	45
Hospital Products	
ROHO® DRY FLOATATION® O.R. Pad	48
Thyroid Pillow	48
Nesting Pad	48
Head/Neck Pillow	48
Interface Dragoure Massurement Customs	
Interface Pressure Measurement Systems	40
XSensor® X3	
XSensor® X3 Display	
XSensor® X3 System	
XSensor® X3 Sensor Pads	
XSensor® X3 LX100	
XSensor® X3 Software	52
Seating Accessories	
Seating Accessories Contour Base	53
The state of the s	
Contour Base	53
Contour Base	53 53
Contour Base Planar Solid Seat Insert ROHO® Heavy Duty Cushion Cover	53 53 54
Contour Base	53 53 54
Contour Base	53 53 54 54
Contour Base Planar Solid Seat Insert ROHO® Heavy Duty Cushion Cover Designer Covers Universal Cushion Cover Special Application Cushions Lumbar Support Cushion	53 53 54 54
Contour Base Planar Solid Seat Insert ROHO® Heavy Duty Cushion Cover Designer Covers Universal Cushion Cover Special Application Cushions Lumbar Support Cushion Sacral Support Cushion	53 54 54 55
Contour Base Planar Solid Seat Insert ROHO® Heavy Duty Cushion Cover Designer Covers Universal Cushion Cover Special Application Cushions Lumbar Support Cushion Sacral Support Cushion MINI-MAX® Cushion	53 54 54 55 55
Contour Base Planar Solid Seat Insert ROHO® Heavy Duty Cushion Cover Designer Covers Universal Cushion Cover Special Application Cushions Lumbar Support Cushion Sacral Support Cushion MINI-MAX® Cushion PACK-IT® Cushion	53 54 54 55 55 55
Contour Base Planar Solid Seat Insert ROHO® Heavy Duty Cushion Cover Designer Covers Universal Cushion Cover Special Application Cushions Lumbar Support Cushion Sacral Support Cushion MINI-MAX® Cushion PACK-IT® Cushion HEAL PAD®	53 54 54 55 55 56
Contour Base Planar Solid Seat Insert ROHO® Heavy Duty Cushion Cover Designer Covers Universal Cushion Cover Special Application Cushions Lumbar Support Cushion Sacral Support Cushion MINI-MAX® Cushion PACK-IT® Cushion	53 54 54 55 55 56 56
Contour Base Planar Solid Seat Insert ROHO® Heavy Duty Cushion Cover Designer Covers Universal Cushion Cover Special Application Cushions Lumbar Support Cushion Sacral Support Cushion MINI-MAX® Cushion PACK-IT® Cushion HEAL PAD® Toilet Seat Cushion	53 54 54 55 55 55 56 56 57
Contour Base Planar Solid Seat Insert ROHO® Heavy Duty Cushion Cover Designer Covers Universal Cushion Cover Special Application Cushions Lumbar Support Cushion Sacral Support Cushion MINI-MAX® Cushion PACK-IT® Cushion HEAL PAD® Toilet Seat Cushion Shower/Commode Cushion ADAPTOR® Pad	53 54 54 55 55 56 56 57
Contour Base Planar Solid Seat Insert ROHO® Heavy Duty Cushion Cover Designer Covers Universal Cushion Cover Special Application Cushions Lumbar Support Cushion Sacral Support Cushion MINI-MAX® Cushion PACK-IT® Cushion HEAL PAD® Toilet Seat Cushion Shower/Commode Cushion	53 54 54 55 55 56 56 57 57

ROHO® DRY FLOATATION® Technology

There are no greater defenses against wounds and disease than those contained in the human body. But often, these defenses need an ally. A partner in the healing process. ROHO DRY FLOATATION technology is that partner.


It creates the ideal healing environment by minimizing soft tissue deformation and maintaining blood flow. And it offers passive compression to improve the efficiency of the metabolic process of the body by reducing edema.

The four principles of ROHO® DRY FLOATATION® Technology




Six Degrees of Freedom

ROHO cushions are constructed of individual cells that move independently, allowing each cell to twist, turn, bend, move up and down, and adapt precisely to the contours of the client.

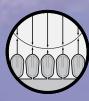

Constant Restoring Forces

When a patient is immersed in a ROHO DRY FLOATATION product, the forces and pressures pushing back are kept equal at all contact points. As the body is immersed, greater contact area is achieved for dispersion of pressure. Thus, the pressure on any one area is minimized.

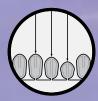

Low Surface Tension

The unique cellular design allows for immersion into the cushion without deforming the tissue, minimizing the chance of skin breakdown.



Low Friction and Shear

As a result of the three previous principles, DRY FLOATATION technology provides a low friction and low shear environment. Friction acts to oppose the direction of motion, or impending motion. Shearing occurs when opposite, but parallel forces meet. The outcome of these forces inhibits blood flow. The slick or non-friction surface of the ROHO cells combined with the independent movement of each cell greatly reduces friction and shear as clients move throughout the day.

How ROHO® DRY FLOATATION® Technology creates a healing environment.

The soft, flexible interconnected air cells that comprise ROHO DRY FLOATATION products are adjustable to provide a customized fit to the client's sitting or lying shape. As the client's condition and body shape changes, the cushion can be adjusted to facilitate blood flow.

WE PROMISE
TO DEVELOP
NEW TECHNOLOGIES
TO IMPROVE
OUR PRODUCTS
AND SERVICES.



"ROHO products have allowed me to be active in my life while staying free from the threat of pressure sores. My ROHO cushion has supported me while winning tennis tournaments in 11 different countries, camping trips with my friends, and even while sitting on our wedding bench during my marriage!"

Michael Foulks, San Diego, California, USA ROHO User Since May 1997



ROHO[®] Seat

The ROHO Group is Shape Fitting Technology. Our technology is unique in its ability to conform to the body's shape. Like fluid, it continuously adjusts with every movement. A person can immerse into the soft cells, and feel almost completely buoyant.

Gels can't do it. Foam can't do it. Even alternative air products cannot do what ROHO technology can do. No other product can more evenly distribute the forces of body weight or reduce the effects of friction or shear upon tissue. ROHO is the optimal choice for preventing or treating pressure sores (ischemic ulcers) and postoperative wounds.

"I've tried all others, none compare to ROHO. Thousands of dollars spent in cushions, and ROHO blows the others away! I ride hard, ROHO rides with me."

Darwin Holmes Modesto, California, USA ROHO User Since July 2003



ROHO® QUADTRO SELECT® Cushion

The QUADTRO SELECT Cushion sets the standard in the industry for overall performance in wheelchair seating. Each and every client will benefit from a customized fit through the simple push of a knob. The revolutionary ISOFLO MEMORY CONTROL® Unit offers shape-fitting capabilities while the client is seated, allowing quick and easy, ondemand adjustment to maximize function. With the built-in stability and simplicity of the QUADTRO SELECT cushion, no longer will you have to sacrifice maximum skin protection to get stability, positioning or convenience.

The QUADTRO SELECT Cushion is available in a 4" (10cm) HIGH PROFILE® model, a 3" (7cm) MID PROFILE™ model and a 2" (6cm) LOW PROFILE® model.

Sizes of the QUADTRO SELECT Cushion come in a wide range of sizes to fit chair widths 12" to 22" / 30cm to 56cm.

Approximate weight of the HIGH PROFILE QUADTRO SELECT is 3.8 lbs. / 1.7 kg, 3.5 lbs. / 1.6 kg for MID PROFILE QUADTRO SELECT and 2.5 lbs. / 1.1 kg for LOW PROFILE QUADTRO SELECT.


36 - Month Limited Warranty.

U.S. Medicare Code: E2624/2625



Patented ISOFLO MEMORY CONTROL® Unit

The SELECT cushions offer versatility with four quadrants controlled by the ISOFLO MEMORY CONTROL Unit.

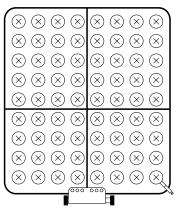

ROHO® CONTOUR SELECT® Cushion

The CONTOUR SELECT Cushion provides the skin protection you expect from a ROHO Group product, with the positioning and stability needed for increased function. The contoured design combined with the simple operation of the ISOFLO MEMORY CONTROL® Unit, provides a new level of postural control and stability that has never been achieved in an air cushion. The contoured design stabilizes the pelvis back in the wheelchair and centers the client comfortably in the middle of the cushion. The ISOFLO MEMORY CONTROL Unit locks air into each of the four quadrants, providing stability by minimizing side-to-side and front-to-back motion.

Sizes of the CONTOUR SELECT Cushion come in a wide range of sizes to fit chair widths 15" to 20" / 38cm to 51cm.

Approximate weight of the CONTOUR SELECT is 3.6 lbs. / 1.6 kg.

36 - Month Limited Warranty

U.S. Medicare Code: E2624/2625





When open, the cushion is inflated through one air valve and the user is adjusted for the lowest bony prominence and positioned for their specific seating needs.

When closed, the ISOFLO MEMORY CONTROL Unit locks the cushion to the required seat position.

Color-coded control helps visually aid in the operation.

ROHO® HIGH PROFILE® Cushion

The ROHO HIGH PROFILE
Cushion is made of 4" (10
cm) interconnected air cells
that provide an adjustable,
low deformation, counter-pressure
environment that assists in the healing of ischemic
ulcers. Appropriate for the high risk client.

*Also available: Dual Compartment can be adjusted independently to increase positioning and stability for either side-to-side or front-to-back control.

The ROHO HIGH PROFILE Cushion comes in a wide range of sizes, to fit chair widths 12" to 28"/ 30cm to 71cm.

Approximate weight of the HIGH PROFILE is 3.8 lbs. / 1.7 kg.

24 - Month Limited Warranty

U.S. Medicare Code: E2622/2623

ROHO® MID PROFILE™ Cushion

The ROHO MID PROFILE
Cushion is made of 3" (7cm)
interconnected air cells that provide an optimal seating environment. The simple adjustment allows the cushions to fit the shape of the client, minimizing tissue deformation and maximizing blood flow. Appropriate for the moderate to high risk client.

The ROHO MID PROFILE Cushion comes in a wide range of sizes, to fit chair widths 13" to 20" / 33cm to 51cm.

Approximate weight of the MID PROFILE is 2.5 lbs / 1.1 kg 24 – Month Limited Warranty

U.S. Medicare Code: E2622/2623

ROHO® LOW PROFILE® Cushion

The ROHO LOW PROFILE
Cushion is made of 2"
(6cm) interconnected air
cells that maximize stability
for the client while protecting from a
potential risk of tissue breakdown. Lightweight
and portable, this cushion is ideal for the active client.

*Also available: LOW PROFILE Dual Compartment can be adjusted independently to increase positioning and stability for either side-to-side or front-to-back control.

The ROHO LOW PROFILE Cushion comes in a wide range of sizes, to fit chair widths 12" to 22" / 30cm to 56cm.

Approximate weight of the LOW PROFILE is 2.5 lbs. / 1.1 kg.

24 - Month Limited Warranty

U.S. Medicare Code: E2622/2623

ROHO® ENHANCER® Cushion

The ENHANCER Cushion uses
AIR IN PLACE® positioning
in a uniquely designed
two-manifold system for
enhanced midline channeling of
the femurs, lateral stability and tissue
protection.

The ENHANCER Cushion comes in a wide range of sizes, to fit chair widths 12" to 20" / 30cm to 51cm.

Approximate weight of the ENHANCER is 3.5 lbs. / 1.6 kg.

24 - Month Limited Warranty

U.S. Medicare Code: E2624/2625

ROHO® Hybrid Elite® Cushion

The ROHO Hybrid Elite Cushion combines the superb skin protection and shape matching ability of ROHO DRY FLOATATION® with the stability of a customized



Jay® contoured foam base. This cushion design allows for improved pelvic alignment and stability for improved position and function while maintaining optimal immersion and envelopment of both the ischial tuberosities and greater trochanters. Available in single and dual valve models.

The ROHO Hybrid Elite Cushion comes in a range of sizes, to fit chair widths 14" to 24" / 35cm to 60cm

Approximate weight of the ROHO Hybrid Elite is 3.5 lbs. / 1.5 kg.

24 - Month Limited Warranty
Medicare Code: E2622/2623

ROHO® nexus SPIRIT® Cushion

The nexus SPIRIT Cushion combines the stability of a contoured foam base with the therapeutic care for which ROHO® DRY FLOATATION® cushions are known. The



contoured foam base is designed to allow for increased stability for transfers and provides exceptional positioning of the pelvis and lower extremities for enhanced sitting posture.

The nexus SPIRIT Cushion comes in a wide range of sizes, to fit chair widths 14" to 20" / 36cm to 51cm.

Approximate weight of the nexus SPIRIT Cushion is 3 lbs. / 1.4 kg.

24 - Month Limited Warranty U.S. Medicare Code: E2622

ROHO® HARMONY® Cushion

By combining a lightweight and comfortable ROHO® AIR FLOTATION™ DEHP-free polyvinylchloride cushion with a pre-contoured foam base, encased within a moisture resistant cover, the client can maximize their independence without the worry of adding unnecessary weight to their mobility system.

The HARMONY Cushion comes in a wide range of sizes, to fit chair widths 14" to 20" / 36cm to 51cm.

Approximate weight of the HARMONY Cushion is 2 lbs. / 1 kg. 250 lbs. / 113.4 kg. weight limit

12 - Month Limited Warranty U.S. Medicare Code: E2607

The ROHO® AirLITE Cushion combines the benefits of

ROHO® AirLITE® Cushion

a sealed air support pad and a durable contoured foam to provide postural support and optimal durability. The non-adjustable pre-set air support pad, uses patented ROHO® AIR FLOATATION technology, while a contoured high density foam base positions the body in a comfortable position, creating one of the most effective and easy-to-use cushions on the market today. And at only 2 pounds, it continues the tradition of lightweight products

The AirLITE cushion fits chair widths 14" to 20" / 35.5 cm to 51 cm.

Approximate weight of the AirLITE cushion is 2.17 lbs. / 1 kg 250 lbs. / 114 kg weight limit

24-Month Limited Warranty
U.S. Medicare Code: E2605

from The ROHO Group.

ROHO® MOSAIC® Cushion

The MOSAIC Cushion is made of 3" (7.5cm) DEHP-free polyvinylchloride interconnected air cells that are easy to inflate



and adjust. Designed to provide a comfortable, functional and stable sitting environment.

The MOSAIC Cushion fits chair widths 16" to 18" / 41cm to 46cm. Available with a standard cover or incontinent cover.

Approximate weight of the MOSAIC Cushion is 0.9 lb. / .4 kg. 250 lbs. / 113.4 kg. weight limit

12 - Month Limited Warranty U.S. Medicare Code: E2601

PostureLITE™ by ROHO®

The PostureLITE™ by ROHO is a self-adjusting wheelchair cushion that combines unparalleled ease of use with comfortable skin protection. It is ideal for anyone who may be at risk of skin breakdown and may have normal sensation.



POSTUREL1CASE POSTUREL2CASE POSTUREL3CASE

PostureLITE comes in three sizes adapt to fit chair sizes ranging from 15" to 20" / 38 cm to 51 cm

POSTUREL1: 16.75 in x 17.25 in (42.5 cm x 44 cm) POSTUREL2: 18.75 in x 17.25 in (47.5 cm x 44 cm) POSTUREL3: 20.00 in x 17.25 in (51.0 cm x 44 cm)

Approximate weight of the PostureLITE Cushion is 1.5 lbs. / 0.6 kg. 350 lbs. / 158 kg. weight limit

12 - Month Limited Warranty U.S. Medicare Code: E2601

LTV Seat® Cushion

You'll Love that Valve!
The valve lets you
customize the seat
to your specific comfort
level. ROHO Shape Fitting
Technology makes the difference. The
LTV Seat Cushion is made of individual,
interconnected air cells that allow air to slowly transfer from chamber to chamber and evenly distribute body
weight pressure. This lightweight Seat Cushion is easy to
fold and can easily be transported from your scooter, to the
car or wherever you may need a cushion for use outside of
a wheelchair.

The LTV Seat Cushion is a comfort cushion. It is not a medical device.

The LTV Seat Cushion is available in:



Approximate weight of the LTV Seat cushion is 1 lb. / .45 kg.

The size of the LTV Seat Cushion is 17" deep x 19" wide / 43cm x 48cm

12 - Month Limited Warranty

Seat Cushion

The following sizes are consistent with most ROHO® DRY FLOATATION® standard adult seat cushions. Please follow these easy steps to create the appropriate model number for one of the DRY FLOATATION seating products.

- Select the appropriate prefix for the product in question. QS, 1R (i.e. Single Comp.), 2R (i.e. Dual Comp.), ENH or CS.
- Select the appropriate number of cells (see below for chart) and place behind the previously selected prefix. (i.e. QS88-This represents a ROHO® QUADTRO SELECT® cushion which is 8 cells wide and 8 cells deep).

Finally, place a C behind the entire letter/ number combination. If the product in question is a ROHO® LOW PROFILE® (either single compartment or dual compartment), place an LPC behind the previous characters.

Example: QS88LPC: LOW PROFILE® QUADTRO SELECT® cushion, 15" x 15" ENH1010C: ENHANCER® cushion, 18" x 18" 1R910LPC: Single Compartment R0H0® LOW PROFILE® cushion, 16" x 18" 2R99C: Dual Compartment R0H0® HIGH PROFILE® cushion, 16" x 16" or 17" x 17"

Sizing Chart

Chair Size (width x depth)	Recommended Number of Cells
12" x 12" / 30 cm x 30 cm	6 wide x 6 deep
12" x 13" / 30 cm x 33 cm	6 wide x 7 deep
13" x 13" / 33 cm x 33 cm	7 wide x 7 deep
14"-15" x 14"-15" / 36-38 cm x 36-38 cm	8 wide x 8 deep
13" x 14"-15" / 33 cm x 36-38 cm	7 wide x 8 deep
13" x 16"-17" / 33 cm x 41-43 cm	7 wide x 9 deep
14"-15" x 16"-17" / 36-38 cm x 41-43 cm	8 wide x 9 deep
14"-15" x 13" / 36-38 cm x 33 cm	8 wide x 7 deep
14"-15" x 14"-15" / 36-38 cm x 36-38 cm	8 wide x 8 deep
14"-15" x 16"-17" / 36-38 cm x 41-43 cm	8 wide x 9 deep
14"-15" x 18"-19" / 36-38 cm x 46-48 cm	8 wide x 10 deep
16"-17" x 14"-15" / 41-43 cm x 36-38 cm	9 wide x 8 deep
16-17" x 16-17" / 41-43 cm x 41-43 cm	9 wide*x 9 deep
16-17" x 18-19" / 41-43 cm x 46-48 cm	9 wide* x 10 deep
16"-17" x 20" / 41-43 cm x 51 cm	9 wide* x 11 deep
18"-19" x 14"-15" / 46-48 cm x 36-38 cm	10 wide* x 8 deep
18"-19" x 16"-17" / 46-48 cm x 41-43 cm	10 wide* x 9 deep
18-19" x 20" / 46-48 cm x 51 cm	10 wide x 11 deep
18"-19" x 18"-19" / 46-48 cm x 48-46 cm	10 wide* x 10 deep
20" x 16"-17" / 51 cm x 41-43 cm	11 wide* x 9 deep
20" x 18"-19" / 51 cm x 46-48 cm	11 wide* x 10 deep
20" x 20" / 51 cm x 51 cm	11 wide* x 11 deep
22" x 18"-19" / 56 cm x 46-48 cm	12 wide* x 10 deep
24" x 18"-19" / 61 cm x 46-48 cm	13 wide x 10 deep
24" x 20" / 61 cm x 51 cm	13 wide x 11 deep
25" x 20" / 64 cm x 51 cm	14 wide x 11 deep
28" x 20" / 69 cm x 51 cm	15 wide x 11 deep

^{*} Wheels with maximum camber or arms with fabric clothing guards may require narrower cushion size for proper fit. Additional sizes may be available - call for a quotation. Not all sizes are available for the QUADTRO SELECT®, CONTOUR SELECT®, ENHANCER® or LTV® Seat Cushion.

For additional information concerning sizing or custom sizes, please contact Customer Care in the U.S. at (800) 851-3449 or outside the U.S. at 1-618-277-9150.

19

Hybrid Elite® Cushion Single Valve Model Standard Sizes

Model	Standard Sizes
1RHY1414C	14" x 14" / 36 cm x 38 cm
1RHY1415C	14" x 15" / 36 cm x 38 cm
1RHY1416C	14" x 16" / 36 cm x 43 cm
1RHY1515C	15" x 15" / 38 cm x 41 cm
1RHY1516C	15" x 16" / 38 cm x 41 cm
1RHY1517C	15" x 17" / 38 cm x 43 cm
1RHY1615C	16" x 15" / 41 cm x 38 cm
1RHY1616C	16" x 16" / 41 cm x 41 cm
1RHY1617C	16" x 17" / 41 cm x 43 cm
1RHY1618C	16" x 18" / 41 cm x 46 cm
1RHY1620C	16" x 20" / 41 cm x 51 cm
1RHY1716C	17" x 16" / 43 cm x 41 cm
1RHY1717C	17" x 17" / 43 cm x 43 cm
1RHY1718C	17" x 18" / 43 cm x 46 cm
1RHY1816C	18" x 16" / 46 cm x 41 cm
1RHY1817C	18" x 17" / 46 cm x 43 cm
1RHY1818C	18" x 18" / 46 cm x 46 cm
1RHY1820C	18" x 20" / 46 cm x 51 cm
1RHY2018C	20" x 18" / 51 cm x 46 cm
1RHY2020C	20" x 20" / 51 cm x 51 cm
1RHY2218C	22" x 18" / 56 cm x 46 cm
1RHY2220C	22" x 20" / 56 cm x 51 cm
1RHY2418C	24" x 18" / 61 cm x 46 cm
1RHY2420C	24" x 20" / 61 cm x 51 cm

Hybrid Elite® Cushion Dual Valve Model Standard Sizes

14" x 14" / 36 cm x 38 cm
14" x 15" / 36 cm x 38 cm
14" x 16" / 36 cm x 43 cm
15" x 15" / 38 cm x 41 cm
15" x 16" / 38 cm x 41 cm
15" x 17" / 38 cm x 43 cm
16" x 15" / 41 cm x 38 cm
16" x 16" / 41 cm x 41 cm
16" x 17" / 41 cm x 43 cm
16" x 18" / 41 cm x 46 cm
16" x 18" / 41 cm x 51 cm
17" x 16" / 43 cm x 41 cm
17" x 17" / 43 cm x 43 cm
17" x 18" / 43 cm x 46 cm
18" x 16" / 46 cm x 41 cm
18" x 17" / 46 cm x 43 cm
18" x 18" / 46 cm x 46 cm
18" x 20" / 46 cm x 51 cm
20" x 18" / 51 cm x 46 cm
20" x 20" / 51 cm x 51 cm
22" x 18" / 56 cm x 46 cm
22" x 20" / 56 cm x 51 cm
24" x 18" / 61 cm x 46 cm
24" x 20" / 61 cm x 51 cm

⁺ May take up to Three weeks for shipment.

nexus SPIRIT® Cushion

Model	Standard Sizes with 2 Support Pad
NS1414C	14" x 14" / 36cm x 36cm
NS1416C	14" x 16" / 36cm x 41cm
NS1515C	15" x 15" / 38cm x 38cm
NS1616C	16" x 16" / 41cm x 41cm
NS1618C	16" x 18" / 41cm x 46cm
NS1717C	17" x 17" / 43cm x 43cm
NS1816C	18" x 16" / 46cm x 41cm
NS1818C	18" x 18" / 46cm x 46cm
NS2018C	20" x 18" / 51cm x 46cm

Note: Sizes listed represent chair size. Actual cushion size is approximately 1.5"/3.8cm deeper to allow for placement between back support canes.

MOSAIC® Cushion

Model (with standard cover)	Model (with incontinent cover)	Fits Chair Size
MOSAIC1616C	MOSAIC1616HD	16" x 16" / 41cm x 41cm
MOSAIC1816C	MOSAIC1816HD	18" x 16" / 46cm x 41cm
		or 16" x 18" / 41cm x 46cm
MOSAIC1818C	MOSAIC1818HD	18" x 18" / 46cm x 46cm

AirLITE® Cushion

Model	Fits Chair Size
AL1414B	14" x 14" / 35.5cm x 35.5cm
AL1515B	15" x 15" / 38cm x 38cm
AL1616B	16" x 16" / 41cm x 41cm
AL1618B	16" x 18" / 41cm x 46cm
AL1716B	17" x 16" / 43cm x 41cm
AL1717B	17" x 17" / 43cm x 43cm
AL1718B	17" x 18" / 43cm x 46cm
AL1720B	17" x 20" / 43cm x 51cm
AL1816B	18" x 16" / 46cm x 41cm
AL1818B	18" x 18" / 46cm x 46cm
AL1820B	18" x 20" / 46cm x 51cm
AL1917B	19" x 17" / 48cm x 43cm
AL1919B	19" x 19" / 48cm x 48cm
AL1920B	19" x 20" / 48cm x 51cm
AL2018B	20" x 18" / 51cm x 46cm

HARMONY® Cushion

Model	Fits Chair Size
H1416C	14" x 16" / 36cm x 41cm
H1616C	16" x 16" / 41cm x 41cm
H1618C	16" x 18" / 41cm x 46cm
H1716C	17" x 16" / 43cm x 41cm
H1718C	17" x 18" / 43cm x 46cm
H1816C	18" x 16" / 46cm x 41cm
H1818C	18" x 18" / 46cm x 46cm
H2018C	20" x 18" / 51cm x 46cm

MINI-MAX® Cushion

Model	Fits Chair Size
MINI1515C	15" x 15" / 38cm x 38cm
MINI1616C	16" x 16" / 41cm x 41cm
MINI1816C	18" x 16" / 46cm x 41cm



You Are An **Original.** So Are We...

There's only one ROHO.

We recognize that alternative air cushions market themselves as our equal. In fact, many use ROHO technology as the benchmark for their own promotion. We accept this and welcome such a large compliment.

We are concerned about the consequences of clients being placed on cushions that do not meet adequate quality and performance levels. It's very flattering to be imitated by other manufacturers, but our biggest concern is the health of our clients. All these cushions look basically the same. What you can't see is the quality and strength of the materials, the engineering, research, technology and testing that goes into a ROHO cushion.

A ROHO cushion is a very refined piece of technology and workmanship. More importantly, its performance has been proven thousands of times over. That is not easily duplicated. There is no guarantee that a cushion that looks like a ROHO is going to perform to the high standards of a ROHO.



ROHO®

Sup

The ROHO Group is proactive. We are advancing. We are learning. We are continuously educating our staff, clinicians and clients. We are always improving. It isn't possible to separate who we are from what we make.

Every ROHO product is the embodiment of our research, testing, education, quality, service, imagination and passionate people. We were founded by innovation, and we will always strive to remain the company that provides the very best in shape fitting solutions.

"I've used ROHO mattress overlays for 10 years. The design gives true pressure relief, prevents pressure ulcers and promotes quick healing. I've used ROHO to heal all levels of pressure ulcers on debilitated patients. It's the only overlay that will give heel pressure relief. The bottom line is, the outcomes are better with the ROHO. I don't worry once a patient is on a ROHO."

Jackie Bangert St. Louis, Missouri, USA Using ROHO Products Since 1995



ROHO® DRY FLOATATION® Mattress Overlay

The ORIGINAL non-powered, adjustable, zoned, pressure-reducing ROHO DRY FLOATATION Mattress Overlay is therapeutic, cost-effective technology at its best. Hundreds of interconnected air cells allow a client with any stage ischemic (pressure) ulcer and/or myocutaneous flan to be



and/or myocutaneous flap to benefit from a custom-adjusted healing environment.

- Friction and shear relief surface
- Zoned. Four sections each 20" x 34" / 51cm x 86cm
- Adjustable. 3.5" (8.9cm) high air cells adjust to the shape and contour of the body
- Non-powered. Requires no external power source
- Lightweight. Weighs eight pounds per section (3.6 kg) or 32 pounds for a complete system (14.4 kg)
- Drainage holes allow fluids to escape away from incontinent client
- Durable neoprene is flame resistant and can be easily cleaned or steam autoclaved
- · Repairable. Extends overall life of mattress
- Covers available

Approximate Weight: 8 lbs. per section / 3.6 kg per section 32 lbs. for a complete unit / 14.4 kg for a complete unit

24 - Month Limited Warranty

U.S. Medicare Code: E0371 Approved for Group II Support Surfaces

ROHO® DRY FLOATATION® Bariatric Kit

Works in conjunction with a standard size
ROHO DRY FLOATATION Mattress Overlay
to widen the mattress for use
with bariatric bed frames.
Add-on sections allow the
bariatric client the therapeutic
benefits of the DRY FLOATATION
Technology. No user weight limit.

Approximate Weight: 4 lbs. / 1.8 kg. (per kit) Standard size: 39" and 48" / Custom sizes available

24 - Month Limited Warranty

U.S. Medicare Code: E0371 Approved for Group II Support Surfaces

ROHO® DRY FLOATATION® Isolette

2.5" (6.5cm) interconnected air cells conform and shape to an infant, providing a low pressure, low shear environment ideal for the at risk infant. This non-powered, adjustable overlay fits standard hospital or home-based isolettes.

Approximate Weight: 6 lbs. / 2.7 kg. 24 - Month Limited Warranty

U.S. Medicare Code: E1399

Reusable Mattress Cover

Waterproof, anti- microbial coated fabric. Offers vapor permeability and can be easily cleaned and/or disinfected. Sized 36" x 80" x 8.5", the reusable cover fits over a



standard single bed mattress, as well as the ROHO® DRY FLOATATION® Mattress System or the PRODIGY Mattress Overlay®. Insert and leveling pads are also available. Please contact Customer Care for details.

6 - Month Limited Warranty

ROHO® PRODIGY® Mattress Overlay® System

Non-powered mattress overlay designed for individuals who are at risk for breakdown and/or have Stage 1 or Stage 2 ischemic ulcers.



- Zoned- Three sections each 27" x 35" / 68cm x 89cm
- Adjustable- 6.5" x 6.5" x 3.5" (16.5cm x 16.5cm x 9cm) air cells adjust to the shape and contour of the body
- · Non-powered- Requires no external power source
- Durable polyvinyl material can be easily cleaned with soap and water

Approximate weight: 3.2 lbs per section (9.6 lbs for a complete unit) / 1.5 kg. per section (4.3 kg for a complete unit) Weight limit: 300 lbs / 136.1 kg

12 - Month Limited Warranty

U.S. Medicare Code: E0197 Approved for Group I Support Surfaces

ROHO® SOFFLEX® 2 Mattress Overlay System

The SOFFLEX 2 is a non-powered, adjustable, and reactive support surface for individuals who require skin/soft tissue protection or have other therapeutic needs. The SOFFLEX 2 provides individualized pressure, shear and friction management through force redistribution to an individual's shape.



- Zoned- Three sections each 36" x 27.5" x 3.5" / 91.5 cm x 70 cm x 9cm
- Adjustable- Air cells adjust to the shape and contour of the body
- Non-powered- Requires no external power source
- Durable flame-resistant polyurethane material can be easily cleaned with soap and water

Approximate weight: 2.25 lbs per section (6.5 lbs for a complete unit) / 1 kg. per section (2.9 kg for a complete unit) Weight limit: 330 lbs / 150 kg

18 - Month Limited Warranty

U.S. Medicare Code: E0371 Approved for Group II Support Surfaces

Sierra™ Gel Overlav

Manufactured for The ROHO Group

- available in the U.S. only



- For the prevention of pressure sores
- One-piece construction requires no assembly
- Easy to clean, low friction surface
- Elastic corner straps to secure overlay to mattress
- Base foundation is a high-density solid one-piece foam constructed core
- Convoluted high-density polyurethane foam topper
- Cover material is Unitek™, a soft nylon urethane coated OFX non-vapor permeable material
- Bladders contain a waterlike, non-toxic gel solution
- Components meet all fire safety and infection control standards

Approximate weight: 46 lbs / 21 kg Weight Limit: 300 lbs / 136 kg

Size: 34" x 78" x 3" / 86.5 cm x 198 cm x 8 cm

One year non-prorated warranty on manufacturer's material and workmanship.

U.S. Medicare Code: E0185 Approved for Group I Support Surfaces

SelectAir® Low Air Loss System

Distributed by The ROHO Group

The SelectAir Low
Air Loss System
is quiet, portable,
lightweight and very
reliable. The SelectAir
mattress and power unit
provide a cost-effective
alternative to other pressure management systems.

- The Firm Mattress Mode helps facilitate nursing care during dressing changes and client repositioning.
- The Upright Patient Mode gives added support when the head of the bed is elevated.
- The Foot Pillow prevents unnecessary skin breakdown due to client migration.

Mattress Configurations: Individual baffles divided into 3 sections; head, seat and foot. Integrated side air bolsters reduce the risk of entrapment.



SelectAir® Max Low Air Loss System

Distributed by The ROHO Group

- available in the U.S. and select international markets.

All the benefits of the *Select*Air plus these special features:

- System Lock prevents accidental adjustments, and automatically resets the system in case of power failure.
- Auto Comfort takes the guesswork out of programming. During set-up, mattress sets optimal pressure according to client's weight distribution.
- Pulsation provides three cycle settings that continually readjust the system to ensure client comfort.

Mattress Configurations: Individual baffles divided into 3 sections; head, seat and foot. Integrated side air bolsters reduce risk of entrapment.

Approximate Weight: Power Unit - 14.3 lbs. / 6.4 kg.

Mattress: 7.6 lbs. / 3.4 kg.

System Weight Limit: 400 lbs. / 181.6 kg.

Zippered Cover (top surface): Waterproof, vapor permeable Custom sizes are available, please contact Customer Support.

One Year Limited Warranty on power unit and soft goods.

U.S. Medicare Code: E0277



Bari Select® Low Air Loss System

Distributed by The ROHO Group

Bariatric clients enjoy ultimate comfort on the Bari Select Mattress. It has two extra inches of depth and a 680-pound weight capacity. It offers more comfort to clients, is less expensive than other fully-integrated systems, and is easier for healthcare workers to operate.

- Quickset allows you to input the client's actual weight, so you can preset the pressure.
- Optional Hand-Held Control allows client to make adjustments as necessary.

Two extra inches of mattress depth, extended widths and specialized software accommodate extra weight capacity.

Mattress Configurations: Individual baffles divided into 3 sections; head, seat and foot. Integrated side air bolsters reduce risk of entrapment.

Approximate Weight: Power Unit - 14.3 lbs. / 6.4 kg.

Mattress: 9.7 lbs. / 4.3 kg.

System Weight Limit: 680 lbs. / 308.7 kg.

Zippered Cover (top surface): Waterproof, vapor permeable

One Year Limited Warranty on power unit, hand control and soft goods. *BariSelect Power Unit is ideal for clients who are between 300 and 680 lbs.

U.S. Medicare Code: E0277



Hybrid Select® with Heal Float® **Low Air Loss System**

Distributed by The ROHO Group - available in the U.S. and Canada

The Hybrid Select with Heal Float is a low air loss mattress system that revolutionizes the overlay market by providing customized



protection. The Hybrid Select with Heal Float combines the power of SelectAir® Low Air Loss Systems with the ultimate skin protection of a ROHO® HIGH PROFILE® Cushion in the heel area to create an optimal environment for skin healing and protection. This system has it all: skin protection, versatility and simplicity!

Size offered: Mattress - 4" height x 80" length x 36" width (10cm height x 213cm length x 91cm width). Cushion - 18" x 16" / 48cm x 41cm. System Weight Limit: 300 lbs. / 136 kg. Warranty: 1 - Year Limited Warranty on Power Unit and soft goods.

2 - Year Limited Warranty on ROHO Cushion. U.S. Medicare Code: E0372 33

Fusion 1K[™] Dynamic Dispersion[™] Mattress Replacement System Distributed by the ROHO Group

The first non-powered therapy mattress that converts to "true" low air loss therapy.

An economical option to the FUSION 1K non-powered therapy.

Dynamic Dispersion™ differs from all other non-powered alternatives. Each air and foam-filled cell communicates to create exceptional envelopment and pressure redistribution.

Dynamic Dispersion

Proven Redistribution Therapy

- Dynamic Dispersion's "closed" technology equalizes and redistributes pressure similar to other proven cellular based technologies.
- The competitor's "open" Self Adjusting Technology design exhausts, therefore the individual foam cells shoulder the majority of the load.
- Since the air is supporting the majority of the load, as opposed to the foam, the system is designed to last years beyond the competition.
- Soft, multi-density foam cradles heels to reduce pressure and shearing in this critical zone.
- High resilient foam, encapsulated in preset air chambers with air flow control. The individual cells communicate to evenly distribute the patient's pressure over the entire therapy zone.
- Top stretch cover is anti-microbial, waterproof and breathable.



Weight Limit: 500 lbs. (227 kg) for system

Warranty: 5-year non-prorated warranty; 1-year warranty on cover


Medicare Code: E0373 Group II pending

Fusion 2K[™] Mattress Replacement System Featuring SelectAir® Low Air Loss Therapy Distributed by The ROHO Group

The first non-powered therapy mattress that converts to "true" low air loss therapy.

Preventative / Therapeutic

 Easily converts from FUSION 1K to 2K in seconds without replacing the mattress.

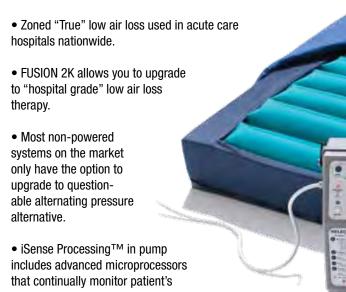

Moisture wicking Patient Moisture Girculates

immersion.



Sizes: 80" x 36" and 84" x 36". 6" height. Approximate Weight: 46 lbs (21 kg) for an average Fusion 2K system

Weight Limit: 300 lbs. (136 kg) for system

Warranty: 5-year non-prorated warranty; 1-year warranty on cover.

Fusion XC[™] Mattress Replacement System Featuring Dynamic Dispersion

Distributed by The ROHO Group


An economical option to the FUSION 1K non-powered therapy.


Dynamic Dispersion™

• Differs from all other non-powered alternatives. Each air and foam-filled cell communicates to create exceptional envelopment and pressure redistribution.

Proven Redistribution Therapy


- The "closed" system envelopes and redistributes pressure similiar to other cellular based technologies.
- The competitor's "open" Self-Adjusting Technology vents making the single foam pieces handle the load individually.
- High resilient foam, encapsulated in preset air chambers with air flow control. The individual cells communicate to evenly distribute the patient's pressure over the entire therapy zone.
- Since the air is supporting the majority of the load, as opposed to the foam, the system is designed to last years beyond the competition.
- Sloped, soft, multi-density foam cradles heels to reduce pressure and shearing in the critical zone.

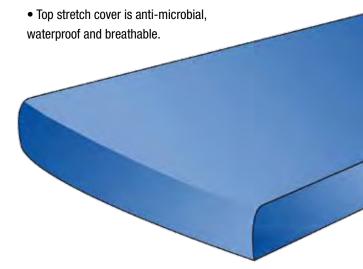

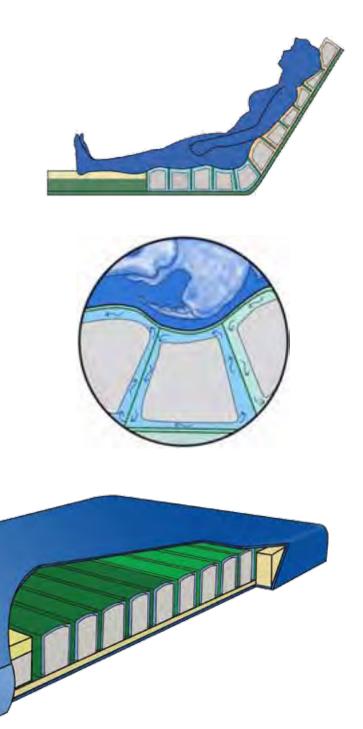

Sizes: 80" x 36" and 84" x 36". 6" height. Weight Limit: 500 lbs. (227 kg) for system

Warranty: 5-year non-prorated warranty; 1-year warranty on cover

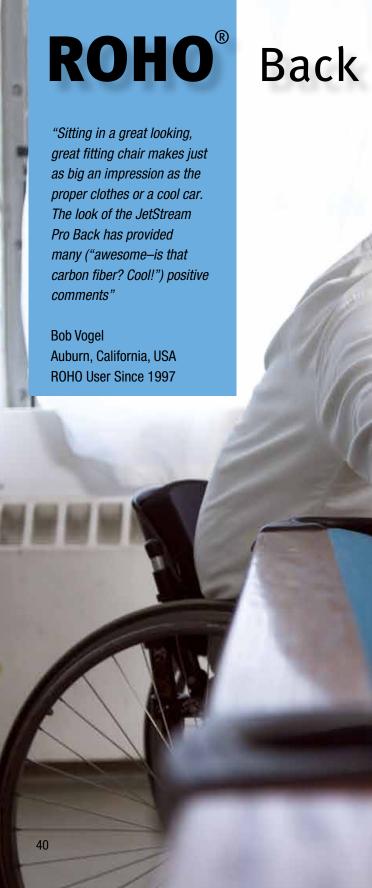

Support Systems



ROHO® AGILITY™ Back System

The New ROHO® AGILITY™ Minimum and Mid Contour Back Systems are the perfect blend of innovative features that give wheelchair users comfortable support without sacrificing function. Innovations that give you peace of mind!

- ROHO AIR FLOATATION™ Insert
- Full Adjustability
- Stability of Design
- Safety You Can Depend On
- Quick Release Hardware
- Easy, Single Hand Removal

One-piece cane clamp fits ¾" to 1" (2 to 3 cm) diameter canes on most wheelchairs and swings around the cane for quick installation.

FLIP-A-BILITY

- Fully adjustable to fit more body shapes.
- The AGILITY™ can be "flipped" to accommodate wider torsos and broader shoulders.
- Full range of Angle, Depth, Width and Vertical adjustment capabilities to address specific positioning needs.



AGILITY™ Back comes standard with:

Back support, mounting hardware, foam pad, lumbar foam support (select models), cover, installation wrench and instruction manual. Adjustable air insert, repair kit and hand pump available on models that include an AIR FLOATATION insert.

AGILITY™ Minimum Contour Back System

Wide range of sizes: Widths 14-20 in. (35.5 cm - 51 cm)

Heights 10-18 in. (25.5 cm - 45.5 cm). Approximate Weight: 2.5lbs - 5.25lbs

Weight Limit: 300 lbs.

Warranty: Back Shell & Hardware - 3 Year Limited Warranty, Air Insert - 18 Month Limited Warranty, Cover & Foam - 6 Month Limited Warranty.

U.S. Medicare Code: E2611, E2613, E2620 and E2615



ROHO® AIR FLOATATION™

Air insert allows individuals to adjust the air volume to their own comfort level and support postural changes.



Laterals also available

AGILITY™ Laterals are designed to provide additional trunk support and positioning for individuals.

The laterals are available in swing-away and fixed brackets with flush or offset mounts.

Swing-away brackets come with curved base with foam pad and cover while fixed brackets come with flat base with foam pad and cover.



Swing Away



Fixed

AGILITY™ Mid Contour Back System

Wide range of sizes: Widths 14-20 in. (35.5 cm - 51 cm) Heights 10-20 in. (25.5 cm - 51 cm).

Approximate Weight: 2.5 lbs - 6 lbs

Weight Limit: 300 lbs.

Warranty: Back Shell & Hardware - 3 Year Limited Warranty, Air Insert - 18 Month Limited Warranty, Cover & Foam - 6 Month Limited Warranty.

U.S. Medicare Code: E2620 and E2615

Back Support Systems

JetStream *Pro*® Back Support System

a ROHO group product Manufactured by ADI

The JetStream *Pro* Back Support System revolutionizes the back support market with a full height back that weighs



only ounces more than standard upholstery. With its lightweight and stable design, this ultralight back is specifically designed to maximize the performance of any mobility base. Sharp, clean lines combined with cutting edge materials make the JetStream *Pro* as appealing to the eye as it is lightweight. The CF (Carbon Fiber) Technology assures the strongest and lightest back available.



JetStream *Pro*® privacy flap standard with purchase



Lateral pads available in either fixed or swingaway models



Adjustable Hardware

Wide range of sizes: Standard Back Height - 16.5" / 42cm,
Mid Back Height - 13" / 33cm, Low Back Height - 10" / 25cm, 14" /
36cm thru 18" / 46cm wide, in 1" / 2.54cm increments
Approximate Weight: 1 lb. / 13.9 oz (.85 kg.) based on 16" / 41cm wide
standard height back, Medium height: 2.25 lb. (1.02 kg.), Approximate
Weight Adjustable Back: 2 lb. (.90 kg.) based on 16" (41cm) wide,
standard height back

24 - Month Limited Warranty

U.S. Medicare Code: E2613 Adjustable & E2611 Standard

Retroback® Back Support System

Ease of use and instant comfort –that's the promise of the new Retroback from The ROHO Group. This unique

new back support system simply clips to the existing sling upholstery to improve comfort and

to the to

support posture. Remove the Retroback by simply lifting it off the upholstery. No tools or hardware required!

The simple yet functional design provides rigidity to compensate for sagging

sling wheelchair back upholstery, while increasing sitting tolerance and reducing lower back fatigue.

Sizes: Models fit wheelchair heights 16" (41cm) to 18" (46cm) and a range of chair widths 15" (38cm) to 20" (51cm). Approximate Weight: Approx. 3 lbs.* (1 kg) based on 18" (46cm) wide back.

*Varies by size.

Weight Limit: 300 lbs.** / 136 kg
**Back support must be properly
sized to client.

12 - Month Limited Warranty
U.S Medicare Code: E2611



Unique Clip On, Clip Off Design

ROHO® Recliner System

Custom designed to the dimensions of the client's chair, this DRY FLOATATION® Recliner System provides back and seating surface protection.

The Recliner System can be manufactured as one single back/seat unit or two separate back/seat cushions.

Approximate average weight of the CUSTOM RECLINER SYSTEM is 12 lbs. / 6 kg.

24 - Month Limited Warranty U.S. Medicare Code:

E2609

ROHO® Hospi

OUR MISSION

To provide superior quality products and support to individuals worldwide who seek skin protection, positioning and comfort solutions.



tal Products



Hospital Products

ROHO® DRY FLOATATION® O.R. Pad

The DRY FLOATATION O.R. Pad is designed for use on operating room tables to promote tissue viability during extended operating procedures. Air cells are 2.5" (6.5 cm) high; drainage holes prevent pooling from caustic solutions.

Approximate Weight for a full system: 42 lbs. / 19 kg. 24 - Month Limited Warranty, Cover not included

ROHO® Thyroid Pillow

The Thyroid Pillow is designed for clients undergoing surgical procedures which may require hyperextension of the neck or other special positions. Able to position client beneath sterile field. Cost-effective due to its r



field. Cost-effective due to its reusability and cleanability.

Approximate Weight: 2 lbs. / 1 kg. **24 - Month Limited Warranty**, Cover not included

ROHO® DRY FLOATATION® Nesting Pad

Designed to simulate the womb environment, the Nesting Pad provides a calmer, better sleeping environment for infants. Multi-shaped air cells create a concave surface and immerse infants into the mattress, conforming to body contours and provides a total contact environment.

Approximate Weight: 7 lbs. / 3 kg. 24 - Month Limited Warranty, Cover not included

ROHO® DRY FLOATATION® Head / Neck Pillow

Provides cushioning and support for the head and neck.

Approximate Weight: 1 lb. / .5 kg. 24 - Month Limited Warranty, Cover not included



Interface Pressure Measurement Systems


XSENSOR® X3

by XSENSOR Technology Corporation - available in the U.S. through The ROHO Group

Interface pressure imaging systems provide objective measurements to assist clinicians in understanding how positioning, support surfaces, and alterations can impact a client's circulation, recovery and well-being. XSENSOR pressure measurement systems are leading the industry in ease of use and portability.

XSENSOR® X3 DISPLAY

by XSENSOR Technology Corporation - available in the U.S. through The ROHO Group

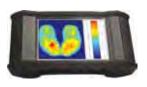



X3 DISPLAY is XSENSOR's answer to the need for improved pressure imaging portability and convenience. This platform is the easiest system to setup, carry, use and share your dynamic pressure images. X3 DISPLAY has pressure imaging software pre-loaded, and autoconfigures to the sensor specifications and calibration files. This platform was designed to meet the needs of clinicians, out-of-lab research and human interface design.


- Handheld portable electronic data logging and data display device
- 6.4 inch 640 x 480 pixel 16 bit resolution LCD screen display with Touchscreen
- . Direct sensor (s) connection to the X3 DISPLAY
- · PC/Laptop USB connection for data downloading and synchronization
- One sensor port connection, which can support up to three individual sensors* at one time
- CompactFlash and Secure Digital memory card ports
- Data synchronization port with two-way communications
- · Connection ports for microphone and headphones

port Internal Lithium Ion batteries provide



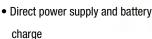

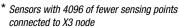

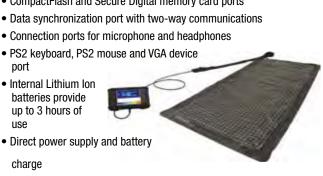



XSENSOR's X3 Technology Series pressure imaging systems offer rehabilitation professionals enhanced plug-and-play functionality and portability. X3 technology provides a dynamic and easy-to-use tool for analyzing and assessing patients. X3 technology provides users with faster processing speeds, increased accuracy and repeatability and enhanced durability.

- · Compact, light-weight electronic data logging device.
- Connects the sensor (s) to a PC/laptop computer USB.
- · Connects to universal power supply.
- The X3 platform is used in a variety of medical, industrial and point-of-sales applications. It provides one data port. Can support up to 3 sensors with the use of an X3 NODE.
- The X3 PRO platform used in multi-sensor or data intensive applications. It provides four data ports and can support up to 12 sensors with the use of X3 NODE.

XSENSOR® X3 SENSOR PADS

by XSENSOR Technology Corporation

- available in the U.S. through The ROHO Group

PX SERIES SENSOR PADS Made with ultra-thin, flexible
sensors, XSENSOR's
systems feature
supple sensor pads
that form to the contours
of uneven surfaces. All sensor
pads are durable and are designed to
withstand rigorous and frequent clinical use, with
slip-resistant polyester covers to ensure easy patient
maneuvering.

- Thin less than 1mm thick
- · Conforming, flexible material conforms to surfaces
- Durable designed to withstand daily patient use
- Holds calibration. May need to be calibrated only 2 to 3 times per year
- High resolution ensures accurate data measurement of smaller body parts

XSENSOR® X3 LX100 High-Accuracy Sensor Series- Now Available

by XSENSOR Technology Corporation

- available in the U.S. through The ROHO Group

LX series sensors are the most advanced and accurate interface pressure sensors offered by XSENSOR.

- · High repeatability
- High level of calibration stability
- In XSENSOR's testing, calibration stability for the LX100 sensor remained within +/- 2% full scale of the initial calibration after five years of accelerated lifetime testing, based on the ASTM 1566-99 (section 9) standard.

XSENSOR® X3 MEDICAL v6.0 Software

by XSENSOR Technology Corporation

- available in the U.S. through The ROHO Group

X3 MEDICAL v6.0 software provides live, dynamic pressure imaging capabilities to clinicians, medical researchers, and durable medical equipment (DME) dealers. Users can choose between two pressure imaging modes. The Seating Assessment Guide uses a simple, four-step process to gather consistent data, compare seating surfaces, and generate detailed reports. The Basic Session mode is a fast way of recording client interface pressures and analyzing data.

XSENSOR Software Version 6.0			
	X3 MEDICAL	X3 DISPLAY	
FEATURES	(installed on PC)	(installed on display)*	
GENERAL			
Capture Static Pressure Image	•	•	
Dynamic Preview Mode	•	•	
Automatic Calibration File Detection	•	•	
Magnifier Window	•	•	
Pressure Range Selection	•	•	
Calibrated or Raw Output	•	•	
Record a Dynamic Pressure Imaging Session	•	•	
Screen Capture	•		
Help Information	•		
Client/Patient Explorer Window	•		
Attach a Subject Image File to the Session	•		
Data Import and Export	•		
Pressure Movie Creation	•		
Session Compare	•		
Sensor Calibration Wizard	Modelled		
Comments/Notes	•	File only, no frame comments	
Video Synch	•		
ANALYSIS AND STATISTICS			
Peak Pressure	•	•	
Contact (Support) Area	•	•	
Average Pressure	•	•	
Mean Reading	•		
Center of Pressure	•		
Sensor Cell Group Analysis	•		
VIEWS			
2D	•	•	
3D	•	•	
Frame compare	•	•	
Histogram	•	•	
Pressure vs. Time	•	•	
Numeric Mode	•		
Configurable Multi-View	•		
* V2 Dieplay packages include V2 N	Indical anthuara a	nd V2 Dianlay anthus	

^{*} X3 Display packages include X3 Medical software and X3 Display software which is pre-loaded from the factory.

Seating Accessories

ROHO® Contour Base

The lightweight, molded base provides a solid foundation for sling seat wheelchairs. Easily fits under the cushion inside the ROHO cushion cover

- · Promotes symmetrical seating
- · Prevents rotation of the femurs
- Provides anti-thrust

Approximate Weight: 1.5 lbs. / 0.7 kg. 12 - Month Limited Warranty



ROHO® Planar Solid Seat Insert

The ROHO Planar Solid Seat Insert provides a solid base for use with ROHO cushioning products. Made of a polyurethane coated birch plywood, the 3/8" (1 cm) insert easily fits under the cushion, inside the ROHO cushion cover. Maximum size - 24" x 24"



- · Eliminates the hammocking effect of sling seats
- · Helps prevent internal rotation of the femurs

Approximate Weight: 2 lbs. / 0.9 kg.

12 - Month Limited Warranty

ROHO® Heavy Duty Cover

Introducing the new ruggedly constructed cover that protects ROHO cushions from rough wear & tear, fluids and the degenerating effects of poor quality ozone. ROHO users can be confident that this new cover will help protect their cushion from rough use and wear and tear. It is also ideal for protection against incontinence with its fluid-resistant material and durable construction. The durable cover design

still allows users to benefit from the therapeutic skin protection of their ROHO cushion.

- · Fluid-resistant, durable material
- Rugged construction
- · Antimicrobial, antibacterial and antifungal
- Double pull zipper allows for easy access to the ISOFLO Memory Control®
- Sizes available to fit most R0H0 cushion models

Available in sizes to fit all standard size ROHO® DRY FLOATATION® seat cushions.


6 - Month Limited Warranty

ROHO® Designer Covers

Personalize your R0H0 $^\circ$ cushion with unique styles that showcase your personality! Custom order to fit your cushion.

91% Polyester / 9% Lycra Top Cover / Snap closure, no zipper


Available for sale at www.therohogroup.com

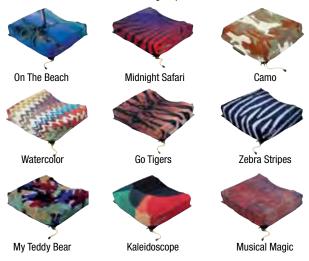

Zebra Stripes

ROHO® Universal Cushion Cover

The cover has a stretch Lycra® top material that is a durable heavy weight fabric, with a soft feel. The side material is composed of a 3-layered knit polyester spacer fabric, which provides superior air circulation that ventilates and cools the user. The fluid proof, solid, non-skid base material offers a

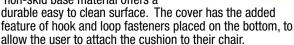

Standard with all HIGH PROFILE®, MID PROFILE™, LOW PROFILE®, QUADTRO SELECT®, MID PROFILE™ QUADTRO SELECT®, LOW PROFILE® QUADTRO SELECT®, CONTOUR SELECT®, ENHANCER®, Custom Recliner System, PACK-IT®, custom and special cushions. Covers can also be custom made.

6 - Month Limited Warranty

^{*}All designer covers are non-returnable. No exchanges.

Special Application Cushions

ROHO® Lumbar Support Cushion

The Lumbar Support cushion is available in a standard model, or can be customized to the needs of the client. Design provides relief from back pain or fatigue. Standard in 4" (10 cm) cell height.

Approximate Weight: 1.25 lbs. / 0.6 kg. **24 - Month Limited Warranty**, Cover not included

ROHO® Sacral Support Cushion

The Sacral Support cushion is custom designed to the client's specification.

Available in 1" (3 cm) to 4" (10 cm) cell heights.

- Designed to help prevent sacral pressure sores
- Encourages proper positioning

Approximate Weight: 1 lb. / 0.5 kg. 24 - Month Limited Warranty, Cover not included

ROHO® MINI-MAX® Cushion

Designed for short term needs of the active client with low risk of skin breakdown, the MINI-MAX cushion provides significant shock and impact absorption, skin protection and stability.



- Minimizes sitting and pressure shearing
- Excellent for sports

Approximate Weight: 2 lbs. / 0.9 kg. **24 - Month Limited Warrantv**

ROHO® PACK-IT® Cushion

Designed for short term needs of active clients, the PACK-IT cushion is excellent for athletic competition, bathing, traveling, or as a scapular, sacral or lumbar support.

- Absorbs shock and impact
- Cell height approx. 1.25"(32 cm)

Approximate Weight: 1 lb. / 0.5 kg. **24 - Month Limited Warranty**

ROHO® HEAL PAD® Cushion

The HEAL PAD offers an adjustable therapeutic environment which aids in the healing of heel and elbow ischemic ulcers.



- Minimizes shear and friction
- · Multiple client use

Approximate Weight: 12 oz. / 96 kg. **24 - Month Limited Warranty**, Cover not included

ROHO® Toilet Seat Cushion

The Toilet Seat cushion is designed for standard size toilet seat and portable bedside commodes. Maximizes stability of the client while protecting from a potential risk of tissue breakdown.



Approximate Weight: 1.25 lbs. / 0.6 kg. 24 - Month Limited Warranty, Cover not included

ROHO® Shower / Commode Cushion

Designed to fit shower / commode chairs, the Shower / Commode Seat cushion offers the therapeutic benefits of ROHO® DRY FLOATATION® technology.

Technology.

Approximate Weight: 1.75 lbs. / 0.8 kg.

24 - Month Limited Warranty, Cover not included

THE ADAPTOR® Pad

The self-contained AIR SPRING® ADAPTOR Pad can be used as an interface on virtually any flat or contoured surface to protect soft tissues.

protect soft tissues.
Especially useful when
cut to fit various surface
shapes such as sides of
wheelchairs, trays, head
and foot rests, inside of
helmets, or anywhere
uncontrolled movement



uncontrolled movement may need protection.

- Absorbs shock
- Protects against shear and friction
- · Enhances soft tissue protection

Approximate Weight: .75 lbs. / 1.3 kg. 24 - Month Limited Warranty, Cover not included

Other replacement parts available:

- Pumps
- Heal Pad Straps
- Inflation Valves
- · Repair Kits
- Cushion Retainer

Special Design Techniques

Custom Variations

Cushions designed to meet a client's size, diagnosis, and needs. All cushions may be custom designed with any combination of the following techniques:

- Vary number of valves / compartments
- Vary cell heights
- Add or delete cell rows beyond dimensions of standard size models
- Remove sections or cells (i.e. cut-outs)
- Ventilation / Drainage openings through base of cushion

Additional product variations:

Lateral Support cushions which attach directly to wheelchair lateral supports.

Custom Molded seating inserts, which fit inside custom molded seats to offer the client increased tissue protection and reduced shock, impact, friction and shear.

Spina Bifida cushions are custom designed based on the needs of the clients with Spina Bifida.

Saddle Cushions, both Western and English models. Custom made to design specifications to minimize shock, impact, friction and shear.

Terms and Policies

Returns: All returns require prior authorization from ROHO, Inc. and are subject to a restocking charge. Before returning your product, contact our Customer Care Department at 1-800-851-3449 toll free in the U.S.A. Outside the U.S.A., contact the nearest ROHO International distributor. See our web site (www.therohogroup.com) for a current list of ROHO International distributors.

Please note that custom ROHO products and XSENSOR products are not returnable. For special return instructions for the Sierra Gel Overlay, contact our Customer Care Department at 1-800-851-3449 toll free in the U.S.A.

Warranty Service Program: Within the U.S., the customer must contact the ROHO, Inc. Warranty Department at 1-800-851-3449 toll free. Note: Any items sent to ROHO, Inc. without prior authorization will not be returned to the customer. Outside the U.S., contact your country's ROHO distributor. For a current list of ROHO distributors, visit our website at www.therohogroup.com, contact ROHO International at 1-618-277-9150 or email us at cc@therohogroup.com.

Be advised that a charge will be assessed for any service provided for a product that does not qualify under the product's warranty.

Note: All seat cushions, except the ADAPTOR® Pad, PostureLITE™, and AirLITE®cushions, include the inflation pump, repair kit and operating instructions.

Notes

ROHO, Inc. has a policy of continual product improvement and reserves the right to amend specifications presented in this document. The current version of this document is available on our website at www.therohogroup.com.

© ROHO, Inc. 2007, 2008, 2009, 2010, 2011, 2012

Following are trademarks and registered trademarks of ROHO, Inc.: ADAPTOR®, AGILITY™, AIR FLOATATION™, AIR IN PLACE®, AIR SPRING®, AIRLITE®, CONTOUR SELECT®, DRY FLOATATION®, ENHANCER®, HARMONY®, HEAL PAD®, HIGH PROFILE®, Hybrid Elite®, ISOFLO MEMORY CONTROL®, JetStream *Pro*®, LOW PROFILE®, MID PROFILE™, MINI-MAX®, MOSAIC®, nexus SPIRIT®, PACK-IT®, PRODIGY Mattress Overlay®, QUADTRO SELECT®, ROHO®, Retroback®, *Select* Protect®, shape fitting technology®, and SOFFLEX®.

PostureLITE™ is a registered trademark of Comfort Concepts, LLC.

BariSelect®, DYNAMIC DISPERSION™, FUSION $1K^m$, FUSION $2K^m$, FUSION XC^m , HealFloat®, HybirdSelect®, iSENSE Processing™, ProtectYou2™ and Select Air® are registed trademarks of MOXI Enterprises, LLC.

XSENSOR® is a registered trademark of XSENSOR Technology Corporation.

Sierra[™] is a trademark of Mason Medical.

Jay® is a trademark of Sunrise Medical.

Ultra Leather™ is a registered trademark of Ultra Fabrics, Inc.

Proven performance, quality & service. When it counts.



100 NORTH FLORIDA AVENUE BELLEVILLE, IL 62221-5429, U.S.A. www.therohogroup.com cc@therohogroup.com U.S. 1.800.851.3449 Fax: 1.888.551.3449

Outside the U.S.: 1.618.277.9150

Fax: 1.618.277.6518

For additional information, please contact: